


BOWIDIC TEST PACK

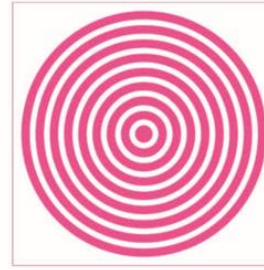

Bowie Dick Test Packs

PRODUCT FEATURES

As mentioned in the European Standard EN 285'a daily steam penetration test should be run in every autoclave. To this end, The Bowie&Dick Autoclave Test Pack is developed. The indicator sheet of this test guarantees colour change from a light blue to an even dark colour and allows for a clear interpretation. The Test Pack complies with ISO 11140-1 (class 2 + 4), EN 867-4 and ISO 11140-4.

The ink is applied to the test sheet trough silk screen printing, which ensures that exactly the same amount of ink is used at all times. The test sheet is laminated to prevent any ink from coming off. The production of all Bowie & Dick Type Autoclave Test Packs is an automatised process under conditioned circumstances. The advantages are that the packs are totally reproducible and the test sheet will stay flat after the sterilization cycle. Full conformity certification for this product has been issued.

| CODE | PRODUCT NAME | FEATURE |
|---------|-------------------|------------|
| ZEO-BW1 | BOWIDIC TEST PACK | 200 SHEETS |
| ZEO-BW1 | BOWIDIC TEST PACK | 400 SHEETS |